

Electric muffle furnaces “N-3 L”, “N-8 L”, “N-13 L”, “N-22 L”, “N-39 L” 1100 °C

FOR TEMPERATURES ADJUSTABLE UP TO 1100 °C.
ELECTRONIC DIGITAL TEMPERATURE CONTROL.
PRECISION ± 2 °C OF THE SET VALUE.
RESOLUTION: 1 DIGIT.




FEATURES

Metal external case. Interior and door made from ceramic fibre, resistant and durable (No asbestos). Heater situated at the side and bottom of the chamber.

CONTROL PANEL

Illuminated mains On/Off switch.
Temperature control with digital display of both the set and actual temperature.
Programmable in steps of 1 °C.
Fitted with a type K probe.

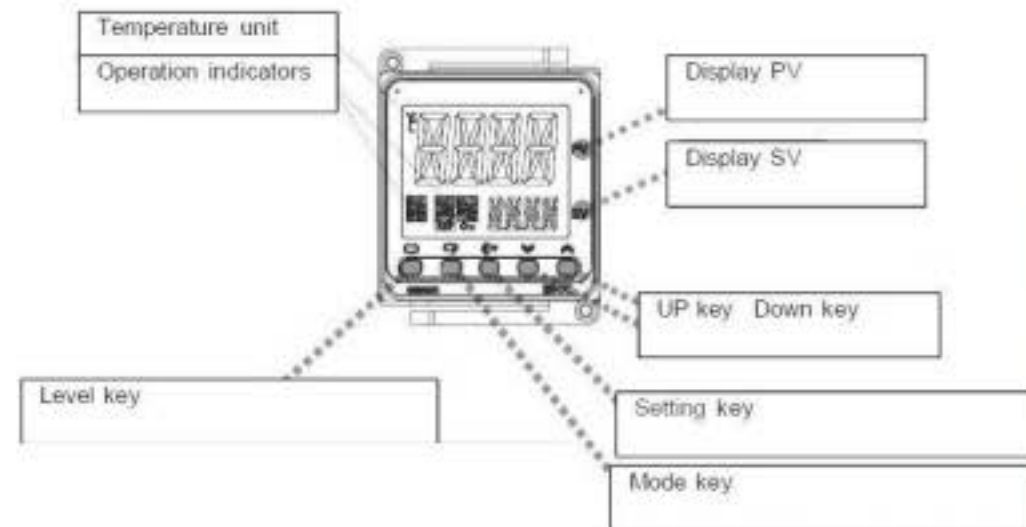

ACCESSORIES

Must be factory installed.

Exterior extractor tube.

Located at the furnace back with a ventilator motor to extract gases and vapours. Gases and Vapours can be extracted outside through the connecting tube. Power: 20 W. Código **2200858**



MODELO	Part No.	Capacity litres	Height / Width / Depht (interior) cm	Height / Width / Depht (exterior) cm	Power W	Voltage V	Weight Kg
N-3 L	2200850	3	11,5 12,5 20 ^o	43 34 47	1700	230	18
N-8 L	2200851	8,2	14 20 30	50 44 53	1800	220	33
N-13 L	2200852	13	18 22,5 36	55 50 70	1800	230	38
N-22 L	2200854	22	15,5 27,5 50	61 60 89	3000	230	58
N-39 L	2200856	39	24 31,5 49,5	74 65 90	6000	400 / 3 N	75

Supplied complete with a refractory ceramic tray as a base and support for material to be assayed.