

TEMPERATURE ADJUSTABLE UP TO 1300 °C FOR MODEL "N-30 L", 1200 °C FOR MODEL "N-40 L" AND 1100 °C FOR MODEL N-80 L.

ELECTRONIC DIGITAL TEMPERATURE CONTROL.

PRECISION ± 2 °C OF THE SET VALUE.

RESOLUTION: 1 DIGIT.

FEATURES

Metal external case with vent at the back of the unit. Interior and door made from ceramic fibre, resistant and durable (No asbestos). Heater situated at the side and bottom of the chamber.



CONTROL PANEL

Illuminated mains On/Off switch.

Temperature control with digital display of both the set and actual temperature.

Programmable in steps of 1 °C.

Fitted with a type K probe.

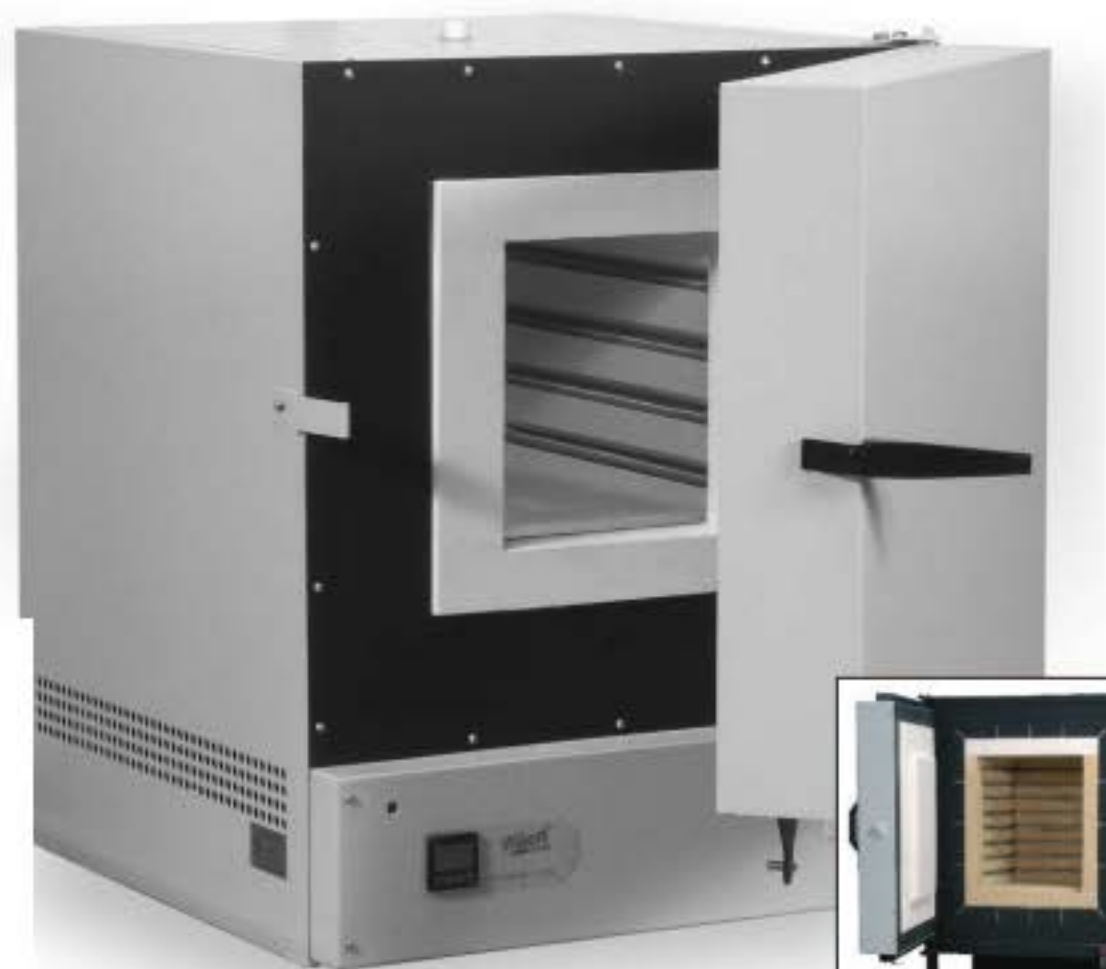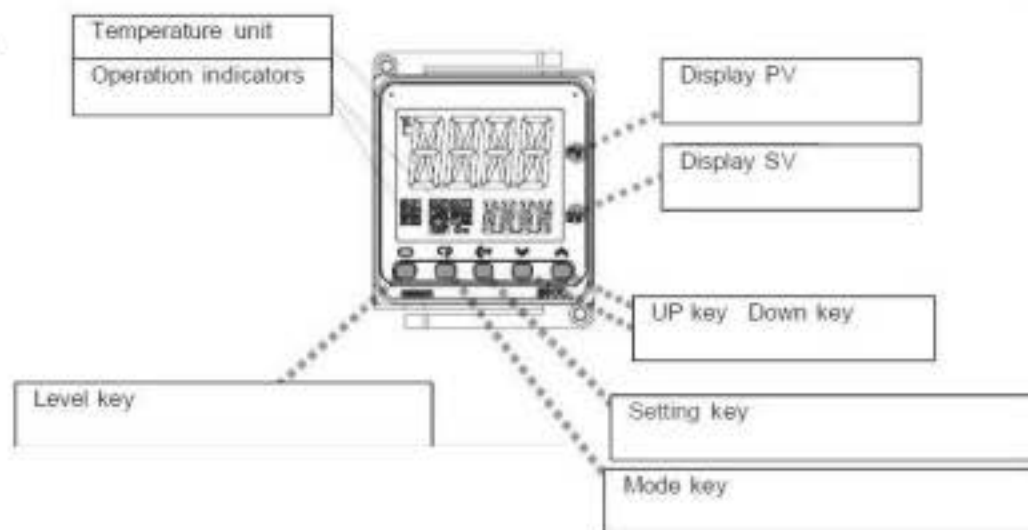

MODEL	Part No.	Capacity litres	Height / Width / Depth (interior) cm	Height / Width / Depth (exterior) cm	Power W	Voltage V	Weight Kg
N-30 L	2200853	30	29 20 44	87 64 84	4600	230	120
N-40 L	2200857	40	29 29 42	87,5 65 83,5	3400	230	105
N-80 L	2200855	80	48 40 40	157 94 98	7500	400 / 3 N	170

Supplied complete with a refractory ceramic tray as a base and support for material to be assayed.